

MITD1 Polyclonal Antibody

Catalog No. E-AB-92084

Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Reactivity	Human,Mouse,Rat
Immunogen	Recombinant fusion protein of human MITD1
Host	Rabbit
Isotype	IgG
Purification	Affinity purification
Conjugation	Unconjugated
Buffer	PBS with 0.01% thiomersal,50% glycerol,pH7.3.

Applications Recommended Dilution


WB	1:500-1:2000
IHC	1:50-1:200

Data

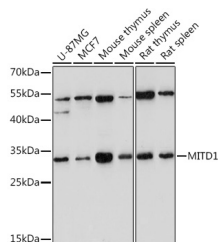

Western blot analysis of extracts of various cell lines using MITD1 Polyclonal Antibody at 1:1000 dilution.

Observed Mw:29kDa
Calculated Mw:29kDa

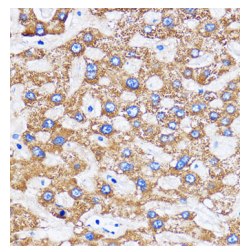

Immunohistochemistry of paraffin-embedded Human liver using MITD1 Polyclonal Antibody at dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

Preparation & Storage

Storage Store at -20°C. Avoid freeze/thaw cycles.

Background

Abscission, the separation of daughter cells at the end of cytokinesis, is effected by endosomal sorting complexes required for transport III (ESCRT-III). The protein encoded by this gene functions as a homodimer, with the N-termini binding to a subset of ESCRT-III subunits and the C-termini binding to membranes. The encoded protein regulates ESCRT-III activity and is required for proper cytokinesis. Several transcript variants encoding different isoforms have been found for this gene.

For Research Use Only